

# d-grip®

## flexible circuit support

The new d-grip® is the perfect airway management support system to prevent load on the airway device from the anaesthesia circuit, monitoring cables, heavy sensors.

### Key features

**Multiple grips** – to fit various sizes of anaesthetic circuit tubing

**Transparent panel** – can be positioned at various angles

**Flexible head section** – allows fine angle and height adjustments

New improved **flexible legs** make adjustment even easier, offering an infinite range of positioning options

d-grip® has been designed to be easy to clean for clinical environments.

**Support channels**  
– for small leads,  
enabling improved  
cable management



Product Code D13001

# d-grip<sup>®</sup> flexible circuit support

Supporting heavy anaesthesia circuits and in-line monitors is essential to avoid upper airway trauma and malpositioning of airway devices leading to potential volatile agent leakage.

The **new d-grip<sup>®</sup>** is the perfect airway management support system to prevent load on the airway device from the anaesthesia circuit, monitoring cables, heavy sensors. It is the perfect complement to **v-gel<sup>®</sup>** and other airway devices.


The transparent panel has various size cut-outs allowing security of a large variety of anaesthetic circuit diameters and even special support channels for small leads such as pulse oximeter cables.

The panel can be removed for cleaning or just wiped over with disinfectant.

The ball and socket joint allows exact positioning of the transparent panel with a simple locking system so provides unlimited positioning options in order to support the circuit at the optimal angle for optimal patient comfort and safety.

The legs are fully flexible to achieve a large scope of height and angle adjustments assisting in finding the perfect positioning required. The legs can even be wound around table edges, drip stands or anaesthetic trolleys for even more fixing options.

**d-grip<sup>®</sup>** gives that extra pair of hands replacing the pile of sandbags and foam wedges and allows veterinary surgeons and nurses to concentrate on treating the patient.

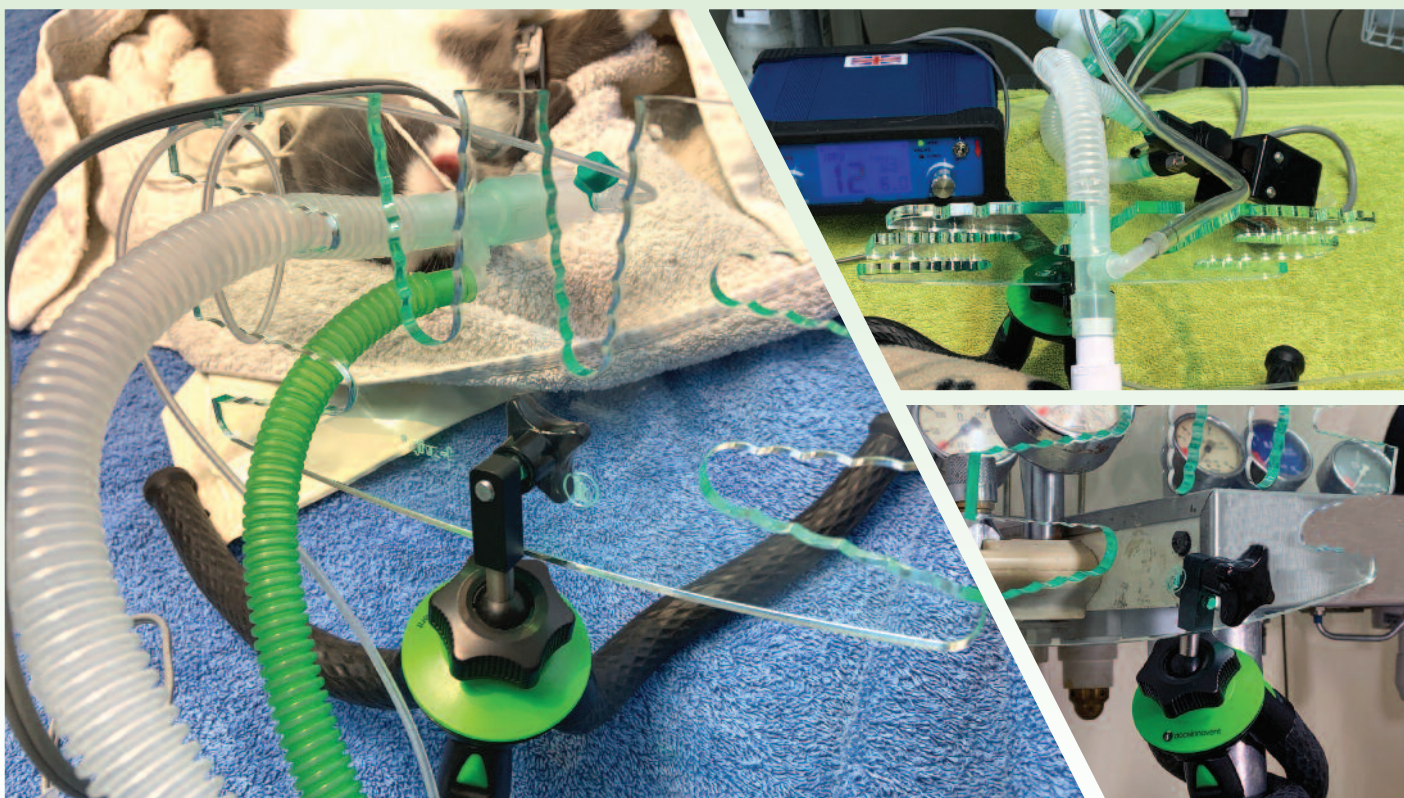

## Advanced veterinary airway management complementary products

The **v-gel<sup>®</sup> Advanced** devices are the first ever true veterinary species specific supraglottic airway devices. **v-gel<sup>®</sup>** devices are available for feline and rabbit species in several sizes. For ordering codes see [www.docsinnovent.com](http://www.docsinnovent.com)



**VetLube<sup>™</sup>** water-based lubricant specifically for veterinary use to assist the insertion of **v-gel<sup>®</sup>** supraglottic airway and other devices. **30ml Spray**  
Order code: **D14002**



To find a distributor in your region please visit [docsinnovent.com/distributors](http://docsinnovent.com/distributors) or contact us directly: +44 (0) 207 649 9071 [info@docsinnovent.com](mailto:info@docsinnovent.com)

